



Manufactured by:

Harvard Medical Devices Ltd. HK

Company: Harvard Medical Devices Ltd. HK

Address: 1002, Railway Plaza, TST, HK



Authorized European Representative:

Company: Share Info Consultant Service LLC Repräsentanzbüro

Address: Heerdter Lohweg 83, 40549 Düsseldorf

Made in China



Kinetik WBP3 UK IB 20220712

wrist blood pressure monitor



- To use the monitor correctly and safely, please read the manual thoroughly.
- Please keep this manual well in order to reference in future.

Instruction Manual

REF WBP3

TABLE OF CONTENTS

INTRODUCTION	2
• General Description	
• Indications for Use	
• Contraindications	
• Measurement Principle	
• Safety Information	
• LCD Display Symbol	
• Monitor Components	
• List	
BEFORE YOU START	13
• Installing and Replacing the Batteries	
• Setting Date, Time, Positioning Symbol and Measurement Unit	
• Select the User ID	
MEASUREMENT.....	22
• Tie the Cuff	
• Start the Measurement	
DATA MANAGEMENT.....	26
• Recall the Records	
• Delete the Records	
INFORMATION FOR USER	30
• Tips for Measurement	
• Maintenance	
ABOUT BLOOD PRESSURE.....	32
• What are systolic pressure and diastolic pressure?	
• What is the standard blood pressure classification?	
• Irregular heartbeat detector	
• Why does my blood pressure fluctuate throughout the day?	
• Why do I get a different blood pressure at home compared to the hospital?	
• Is the result the same if measuring on the right wrist?	
TROUBLESHOOTING	35
SPECIFICATIONS	37
MANUFACTURER INFORMATION.....	39
EMC GUIDANCE	40

♥ General Description

Thank you for selecting Kinetik wrist type blood pressure monitor (WBP3). The monitor features blood pressure measurement, pulse rate measurement and the result storage. The design provides you with two years of reliable service.

Readings taken by the WBP3 are equivalent to those obtained by a trained observer using the cuff and stethoscope auscultation method.

This manual contains important safety and care information, and provides step by step instructions for using the product.

Read the manual thoroughly before using the product.

Features:

- Systolic blood pressure
- Diastolic blood pressure
- Pulse rate
- 60 records for each user

♥ Indications for Use

The Kinetik Wrist Blood Pressure Monitor is a digital monitor intended for use in measuring blood pressure and heartbeat rate with a wrist circumference ranging from 13.5cm to 21.5 cm (about 5 $\frac{1}{3}$ "-8 $\frac{1}{2}$ ").

It is intended for adult indoor use only.

♥ Contraindications

- 1.The device should not be used by any person who may be suspected of,or is pregnant .
- 2.The device is not suitable for use on patients with implanted,electrical devices, such as cardiac pacemakers, defibrillators.

♥ Measurement Principle

This product uses the Oscillometric Measuring Method to detect blood pressure. Before every measurement, the unit establishes a "zero point" equivalent to the atmospheric pressure. Then it starts inflating the cuff. Meanwhile, the unit detects pressure oscillation generated by beat-to-beat pulsatile, which is used to determine the systolic pressure and diastolic pressure as well as pulse rate.

♥ **Safety Information**

The signs below might be in the user manual, labeling or other component. They are the requirement of standard and using.

	Read the instructions (actual symbol colours are white on a blue background).		Symbol for "TYPE BF APPLIED PARTS"
CE0197	Symbol for "COMPLIES WITH MDD 93/42/EEC REQUIREMENTS"		Symbol for "ENVIRONMENT PROTECTION - Electrical waste products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice"
	Symbol for "MANUFACTURER"		
SN	Symbol for "SERIAL NUMBER"		Symbol for "DIRECT CURRENT"
	Symbol for "RECYCLE"		Symbol for "Authorised Representative in the European Community"
	Symbol for "MANUFACTURE DATE"		Caution: These notes must be observed to prevent any damage to the device.

⚠ **CAUTION**

- * This device is intended for adult use in homes only.
- * The device is not suitable for use on neonatal patients, pregnant women, patients with implanted, electrical devices, patients with pre-eclampsia, premature ventricular beats, atrial fibrillation, peripheral, arterial disease and patients undergoing intravascular therapy or arterio-venous shunt or people who received a mastectomy. Please consult your doctor prior to using the unit if you suffer from illnesses.
- * The device is not suitable for measuring the blood pressure of children. Ask your doctor before using it on older children.
- * The device is not intended for patient transport outside a healthcare facility.
- * The device is not intended for public use.
- * This device is intended for no-invasive measuring and monitoring of arterial blood pressure. It is not intended for use on extremities other than the wrist or for functions other than obtaining a blood pressure measurement.
- * Do not confuse self-monitoring with self-diagnosis. This unit allows you to monitor your blood pressure. Do not begin or end medical treatment without asking a physician for treatment advice.
- * If you are taking medication, consult your physician to determine the most appropriate time to measure your blood pressure. Never change a prescribed medication without consulting your physician.
- * Do not take any therapeutic measures on the basis of a self measurement. Never alter the dose of a medicine prescribed by a doctor. Consult your doctor if you have any question about your blood pressure.
- * When the device was used to measure patients who have common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation, the best result may occur with deviation. Please consult your physician about the result.
- * Warning: Do not apply the cuff over a wound; otherwise it can cause further injury.
- * Do not inflate the cuff on the same limb which other monitoring ME equipment is applied around simultaneously, because this could cause temporary loss of function of those simultaneously-used monitoring ME equipment.

INTRODUCTION

CAUTION

- * On the rare occasion of a fault causing the cuff to remain fully inflated during measurement, open the cuff immediately. Prolonged high pressure (cuff pressure > 300mmHg or constant pressure > 15mmHg for more than 3 minutes) applied to the wrist may lead to an ecchymosis.
- * Please check that operation of the device does not result in prolonged impairment of patient blood circulation.
- * When measurement, please avoid compression or restriction of the connection tubing.
- * The device cannot be used with HF surgical equipment at the same time.
- * The ACCOMPANYING DOCUMENT shall disclose that the SPHYGMOMANOMETER was clinically investigated according to the requirements of ISO 81060-2:2013.
- * To verify the calibration of the AUTOMATED SPHYGMOMANOMETER, please contact the manufacturer.
- * This device is contraindicated for any female who may be suspected of, or is pregnant. Besides providing inaccurate readings, the effects of this device on the fetus are unknown.
- * Too frequent and consecutive measurements could cause disturbances in blood circulation and injuries.
- * This unit is not suitable for continuous monitoring during medical emergencies or operations. Otherwise, the patient's wrist and fingers will become anaesthetic, swollen and even purple due to a lack of blood.
- * When not in use, store the device in a dry room and protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on the storage case.
- * This device may be used only for the purpose described in this booklet. The manufacturer cannot be held liable for damage caused by incorrect application.
- * This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in this booklet.
- * The maximum temperature that the applied part can be achieved is 42.5 °C while the environmental temperature is 40 °C.
- * The equipment is not AP/APG equipment and not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
- * Warning: No servicing/maintenance while the ME equipment is in use.
- * The patient is an intended operator.

INTRODUCTION

CAUTION

- * The patient can measure data and change batteries under normal circumstances and maintain the device and its accessories according to the user manual.
- * To avoid measurement errors, please avoid the condition of strong electromagnetic field radiated interference signal or electrical fast transient/burst signal.
- * The blood pressure monitor and the cuff are suitable for use within the patient environment. If you are allergic to polyester, nylon or plastic, please don't use this device.
- * During use, the patient will be in contact with the cuff. The materials of the cuff have been tested and found to comply with requirements of ISO 10993-5:2009 and ISO 10993-10:2010. It will not cause any potential sensitization or irritation reaction.
- * If you experience discomfort during a measurement, such as pain in the wrist or other complaints, press the START/STOP button to release the air immediately from the cuff. Loosen the cuff and remove it from your wrist.
- * If the cuff pressure reaches 40 kPa (300 mmHg), the unit will automatically deflate. Should the cuff not deflate when pressures reaches 40 kPa (300 mmHg), detach the cuff from the wrist and press the START/STOP button to stop inflation.
- * Before use, make sure the device functions safely and is in proper working condition. Check the device, do not use the device if it is damaged in any way. The continuous use of a damaged unit may cause injury, improper results, or serious danger.
- * Do not wash the cuff in a washing machine or dishwasher!
- * The service life of the cuff may vary by the frequency of washing, skin condition, and storage state. The typical service life is 10000 times.
- * The service life of this product is 2 years.
- * Manufacturer will make available on request circuit diagrams, component part lists, descriptions, calibration instructions, etc., to assist to service personnel in parts repair.
- * The operator shall not touch output of batteries and the patient simultaneously.

INTRODUCTION

CAUTION

- * Cleaning :Dust environment may affect the performance of the unit. Please use the soft cloth to clean the whole unit before and after use. Don't use any abrasive or volatile cleaners.
- * The device doesn't need to be calibrated within two years of reliable service.
- * If you have any problems with this device, such as setting up, maintaining or using, please contact the SERVICE PERSONNEL of Harvard Medical Devices. Don't open or repair the device by yourself in the event of malfunctions. The device must only be serviced, repaired and opened by individuals at authorized sales/service centers.
- * Please report to Harvard Medical Devices if any unexpected operation or events occur.
- * Keep the unit out of reach of infants, young children or pets to avoid inhalation or swallowing of small parts. It is dangerous or even fatal.
- * Be careful to strangulation due to cables and hoses, particularly due to excessive length.
- * At least 30 min required for ME equipment to warm from the minimum storage temperature between uses until it is ready for intended use. At least 30 min required for ME equipment to cool from the maximum storage temperature between uses until it is ready for intended use.
- * This equipment needs to be installed and put into service in accordance with the information provided in the ACCOMPANYING DOCUMENTS;
- * Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies can affect this equipment and should be kept at least a distance d away from the equipment. The distance d is calculated by the MANUFACTURER from the 80MHz to 5.8 GHz column of Table 4 and Table 9 of IEC 60601-1-2:2014, as appropriate.
- * Please use ACCESSORIES and detachable partes specified/ authorised by MANUFACTURE. Otherwise, it may cause damage to the unit or danger to the user/patients.
- * There is no luer lock connectors are used in the construction of tubing, there is a possibility that they might be inadvertently connected to intravascular fluid systems, allowing air to be pumped into a blood vessel.
- * Please use the device under the environment which was provided in the user manual. Otherwise, the performance and lifetime of the device will be impacted and reduced.

INTRODUCTION

CAUTION

- * When using this device, please pay attention to the following situation which may interrupt blood flow and influence blood circulation of the patient, thus cause harmful injury to the patient: connection tubing kinking too frequent and consecutive multiple measurements; the application of the cuff and its pressurization on any wrist where intravascular access or therapy, or an arterio-venous (A-V) shunt, is present; inflating the cuff on the side of a mastectomy.
- * Please dispose of ACCESSORIES, detachable parts, and the ME EQUIPMENT according to the local guidelines.


♥ LCD Display Symbol



SYMBOL	DESCRIPTION	EXPLANATION
SYS	Systolic pressure	High blood pressure
DIA	Diastolic pressure	Low blood pressure
PULSE /min	Pulse display	Pulse in beats per minute
	Motion indicator	Motion may result in an inaccurate measurement
 +Lo	Low battery	Batteries are low and need to be replaced


SYMBOL	DESCRIPTION	EXPLANATION
kPa	kPa	Measurement unit of the blood pressure (1kPa=7.5mmHg)
mmHg	mmHg	Measurement unit the blood pressure (1mmHg=0.133kPa)
♥	Irregular heartbeat	Blood pressure monitor is detecting an irregular heartbeat during measurement.
00:00	Current time	Month/Day,Hour/Minute
AVG	Average value	The average value of latest 3 groups blood pressure value.
♥	Heartbeat	Blood pressure monitor is detecting a heartbeat during measurement.
Q	Memory	Indicate it is in the memory mode
1	User 1	Start measurement and save the results for user 1
2	User 2	Start measurement and save the results for user 2
▬	Blood pressure level indicator	Indicate the blood pressure level
👤	Positioning symbol	The correct position gives you an accurate measurement.

♥ Monitor Components

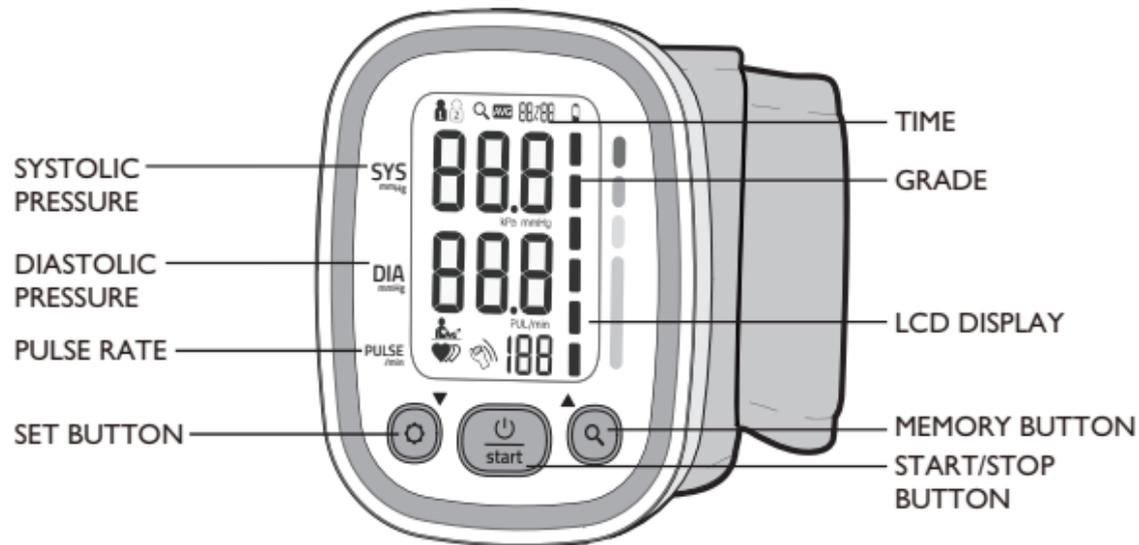

Component list pressure measuring system:

1. PCBA;
2. Air Pipe;
3. Pump;
4. Valve;
5. Cuff.

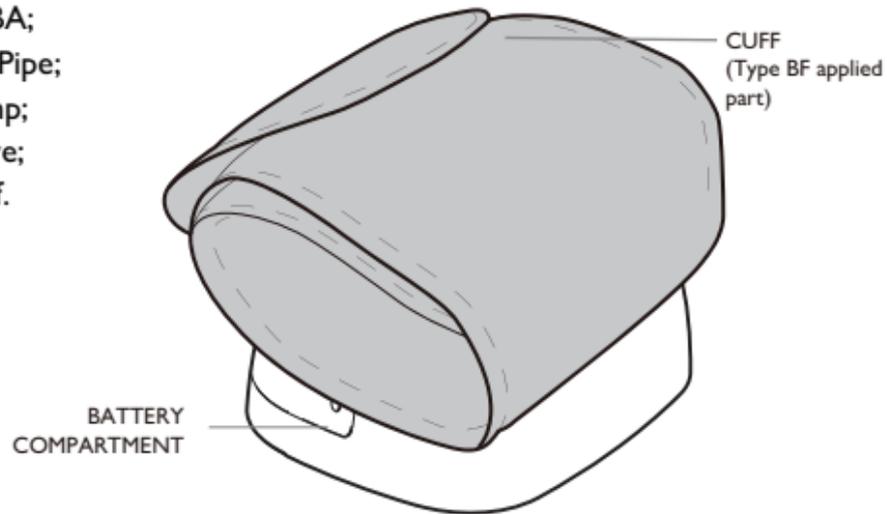


♥ List

- 1) Blood Pressure Monitor WBP3
- 2) 2×AAA batteries
- 3) User manual
- 4) Quick start guide

♥ Installing and Replacing the Batteries

- Slide off the battery cover.
- Install the batteries by matching the correct polarity, as shown below. Always use the correct battery type (2×AAA batteries.)
- Replace the cover.

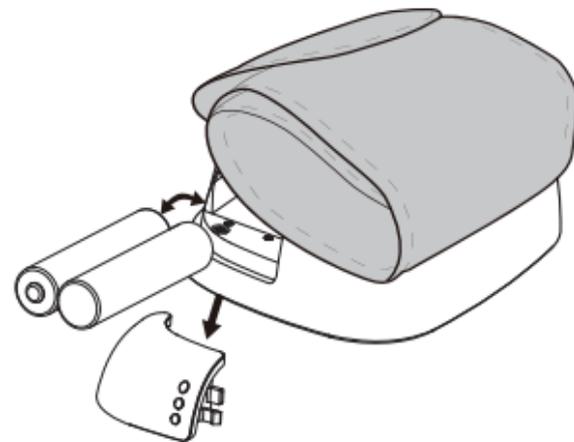

Replace the batteries whenever the below happen

- The  shows
- The display is dim.
- The display does not light up


CAUTION

Do not use new and used batteries together.
 Do not use different types of batteries together.
 Do not dispose the batteries in fire. Batteries may explode or leak.
 Remove batteries if the device is not likely to be used for some time.
 Worn batteries are harmful to the environment. Do not dispose with daily garbage.
 Remove the old batteries from the device following your local recycling guidelines.

♥ Setting Date, Time, Positioning Symbol and Measurement Unit

It is important to set the clock before using your blood pressure monitor, so that a time stamp can be assigned to each record that is stored in the memory. (The setting range of year is 2016 ~ 2056. Time format: 24H)

1. When the monitor is off, press "SET" button, it will display the time. Then press and hold "SET" button to enter the mode for year setting.

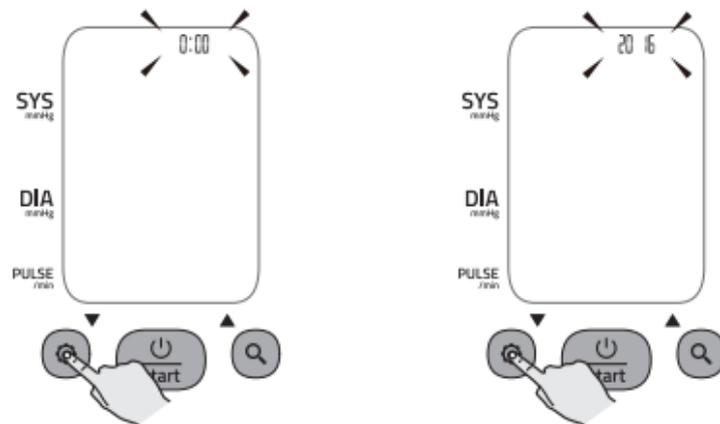


Notes:


1. During the process of setting, you can press “START/STOP” button to stop setting at any time.
2. If there is no operation during the process of setting, it will turn off.

2. Press “MEM” button to change the [YEAR]. Each press will increase the numeral by one in a cycling manner.


3. When you get the right year, press “SET” button to confirm and it will turn to next step.

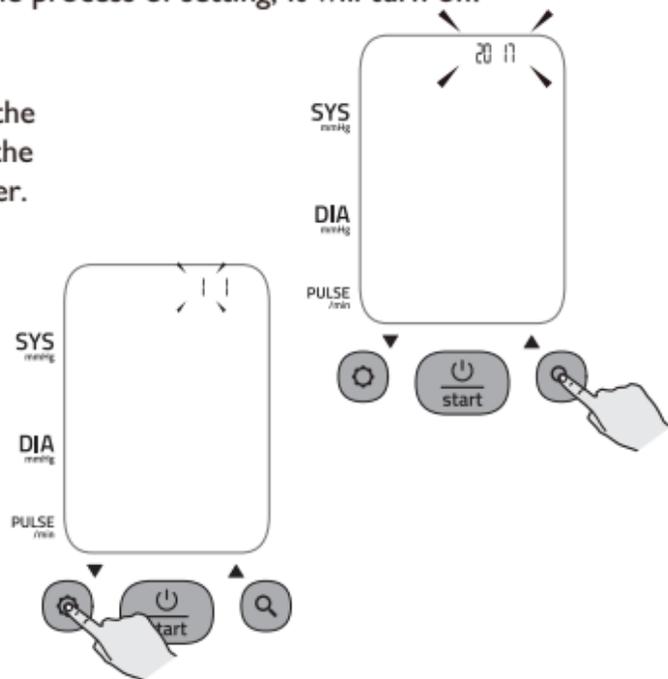

4. Repeat step 2 and 3 to confirm [MONTH] and [DAY].

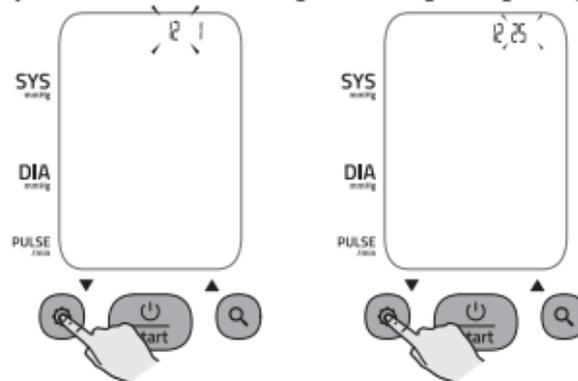


5. Repeat step 2 and 3 to confirm [HOUR] and [MINUTE].

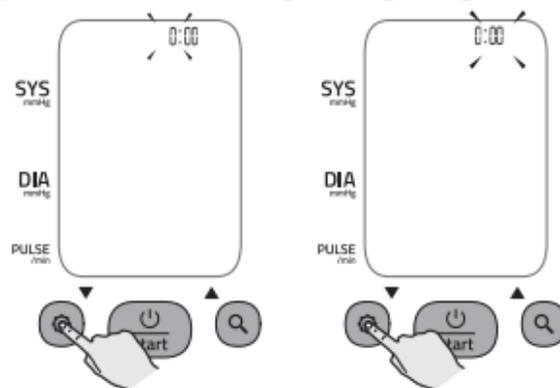


BEFORE YOU START

- Repeat step 2 and 3 to confirm [Positioning symbol] on or off.

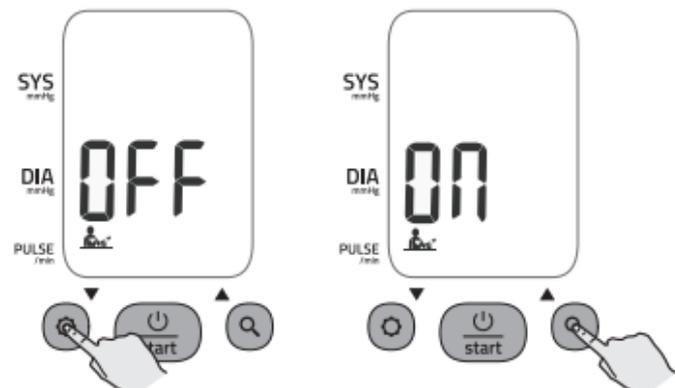


- Repeat step 2 and 3 to confirm the measurement unit.

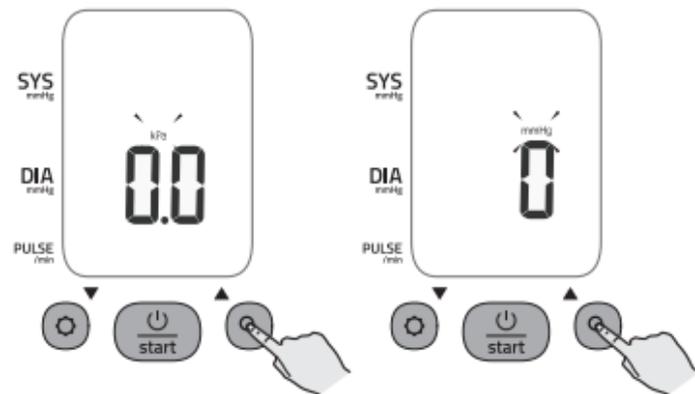

BEFORE YOU START


- After confirming the measurement unit, the LCD will display "do nE", then display all the settings you have done one by one and then the monitor will turn off.

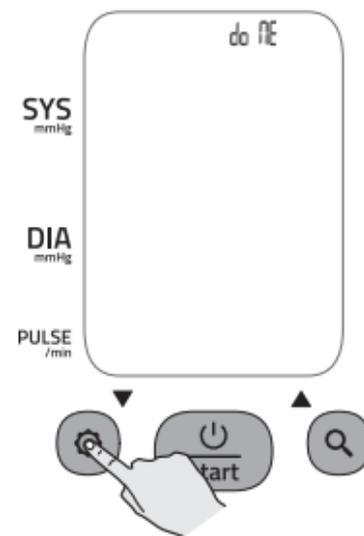


♥ Select the User ID

Before you start the measurement, please select the desired user ID first.

1. When the blood pressure monitor is off, press and hold “MEM” button until the user ID blinks. Then press “MEM” button to change the user ID between user 1 and user 2. Press “SET” button to confirm your selection.

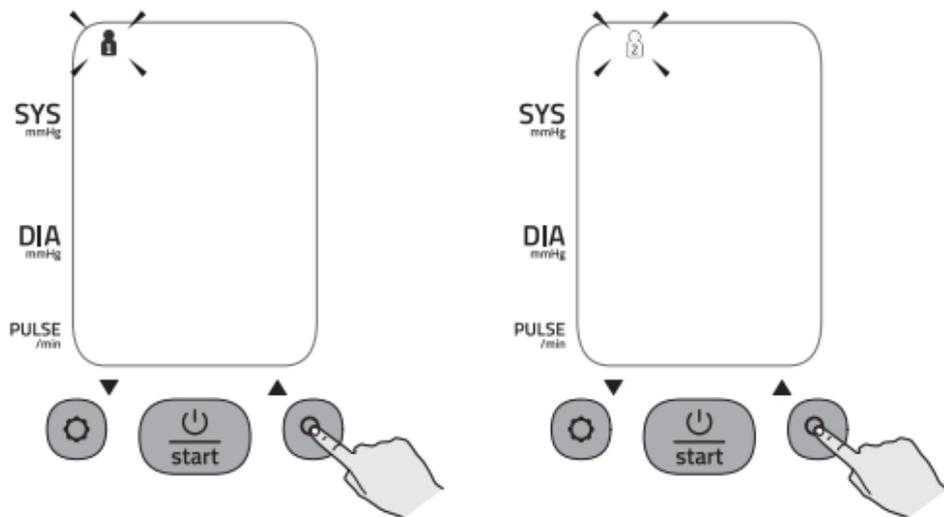

2. After confirming the user ID, the LCD will display “User ID+donE” and then turn off. Then you can start your measurement now.

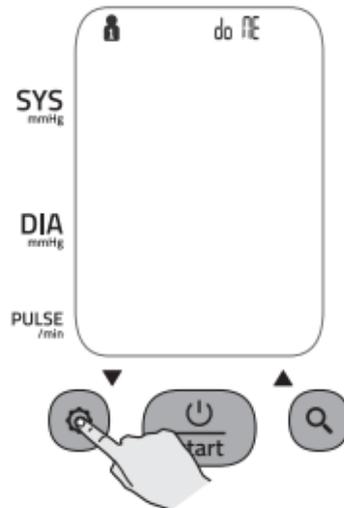

♥ Tie the Cuff

1. Remove all accessories (watch, bracelet, etc) from your wrist. If your physician has diagnosed you with poor circulation in your wrist, use the other one.
 2. Roll or push up your sleeve to expose the skin.
 3. Apply the cuff to your wrist with your palm facing up.
 4. Position the edge of the cuff about 1cm~1.5cm from wrist joints.
 5. Fasten the wrist cuff around your wrist, leaving no extra room between the cuff and your skin. If the cuff is too loose, the measurement will not be accurate.
 6. Sit comfortably with your tested wrist resting on a flat surface. Place your elbow on a table so that the cuff is at the same level as your heart. Turn your palm upwards. Sit upright in a chair, and take 5-6 deep breaths.
 7. Patients with Hypertension:
The middle of the cuff should be at the level of the right atrium of the heart; Before starting measurement, please sit comfortably with legs uncrossed, feet flat on the floor, back and wrist supported.
- Rest for 5 minutes before measuring.
 - Wait at least 3 minutes between measurements. This allows your blood circulation to recover.
 - Take the measurement in a silent room.
 - The patient must relax as much as possible and do not move and talk during the measurement procedure.
 - The cuff should maintain at the same level as the right atrium of the heart.
 - Do not cross your legs and keep your feet on the ground.
 - Keep your back against the backrest of the chair.
 - For a meaningful comparison, try to measure under similar conditions. For example, take daily measurements at approximately the same time, on the same wrist, or as directed by a physician.




♥ Start the Measurement

1. When the monitor is off, press "START/STOP" button to turn on the monitor, and it will finish the whole measurement. (Take user 1 for example.)




LCD display

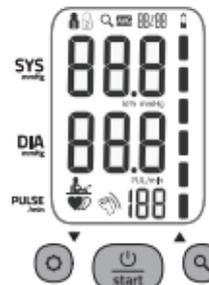

Inflating and measuring.



Adjust the zero.

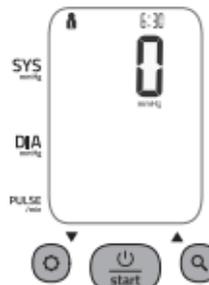

Display and save the result. The year, date and time will display alternately.



MEASUREMENT

2. Press “START/STOP” button to power off, otherwise it will turn off within 1 minute.




Tips:


1. There are two users in total. Each user has 60 records.
2. You can press “START/STOP” button to stop measuring during the process of the measurement at any time.
3. If the measurement result is out of the measurement range (SYS: 60mmHg to 230mmHg; or DIA: 40mmHg to 130mmHg; or Pulse: 40-199 pulse/minute), the LCD will display “out”.

Notes: If you have set the [Positioning Symbol] on, when you start the measurement, the blood pressure monitor will detect the position first. The wrist must beat the angle between 30° and 45°. If it's out of this angle, the LCD won't start any measurements and will display  + Err until you have the correct position.

MEASUREMENT

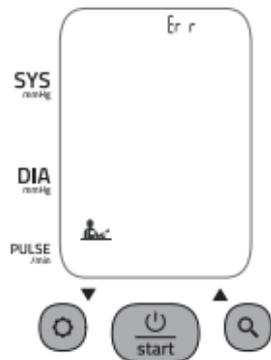

If you get the correct position, it will display  for about 3 seconds, and then finish the whole measurement.

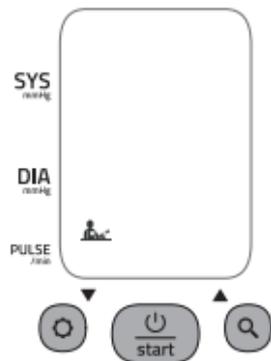

♥ Recall the Records

1. When the monitor is off, press “ MEM ” button to show the average value of the latest three measurement records. If the records are less than 3 groups, it will display the latest record instead.
2. Press “ MEM ” button or “ SET ” button to get the record you want.



The record number, date and time of the measurement records will display alternately.

3. If you want to check the other user's measurement records, please press “ START/STOP ” button to turn off the blood pressure monitor. Then press and hold “ MEM ” button to enter the selecting user mode, press “ MEM ” again to change the user, when the desired user ID blinks, press “ SET ” button to confirm. Then press “ MEM ” button to check the records of the selected user.



CAUTION


The most recent record (1) is shown first. Each new measurement is assigned to the first (1) record. All other records are pushed back one digit (e.g., 2 becomes 3, and so on), and the last record (60) is dropped from the list.

♥ Delete the Records

If you did not get the correct measurement, you can delete all the results for the selected user by following steps.


(Take user 1 for example.)

1. In the memory mode, hold pressing "MEM" button about 3 seconds, the flash display "dEL All" + User ID will show.

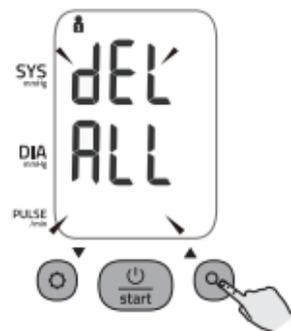


2. Press "SET" button to confirm deleting, the LCD displays "User ID+dEL donE" and then turns off.

Note: To exit out of delete mode without deleting any records, press "START/STOP" button before pressing "SET" button to confirm any delete commands.

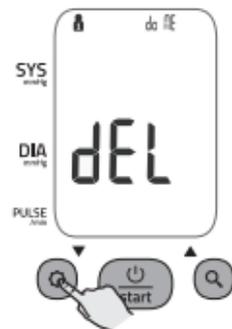

3. If there is no record, when you press "MEM" button to check the record, the display will show.

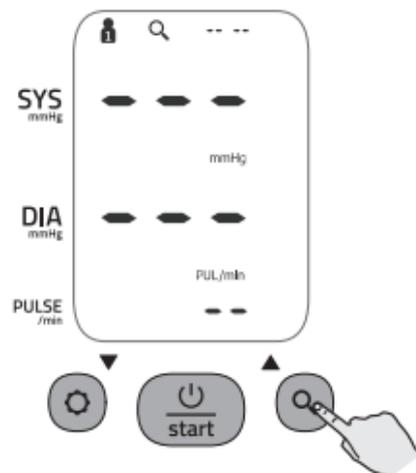

♥ Tips for Measurement

Measurements may be inaccurate if taken in the following circumstances.



Within 1 hour
after dinner or drinking



Immediate measurement
after tea, coffee, smoking



Within 20 minutes
after taking a bath



When talking or
moving your fingers



In a very cold environment



When you want to
discharge urine



♥ Maintenance


To obtain the best performance, please follow the instructions below.



Put in a dry place and
avoid the sunshine



Avoid touching water,
clean it with a dry cloth in case.

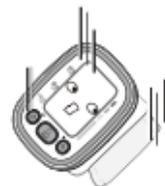

Avoid intense shaking and collision.



Avoid dusty and unstable-
temperature environment.




Using wet cloths to remove dirt.



Avoid washing the cuff

♥ **What are systolic pressure and diastolic pressure?**

When ventricles contract and pump blood out of the heart, blood pressure reaches its maximum value, the highest pressure in the cycle is known as systolic pressure. When the heart relaxes between heartbeats, the lowest blood pressure is diastolic pressure.

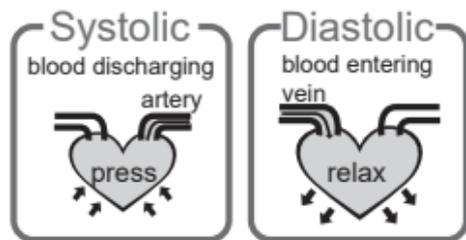


♥ **What is the standard blood pressure classification?**

The blood pressure classification published by World Health Organization(WHO) and International Society of Hypertension (ISH) in 1999 is as follows:

CAUTION

Only a physician can tell your normal BP range. Please contact a physician if your measuring result falls out of the range. Please note that only a physician can tell whether your blood pressure value has reached a dangerous point.

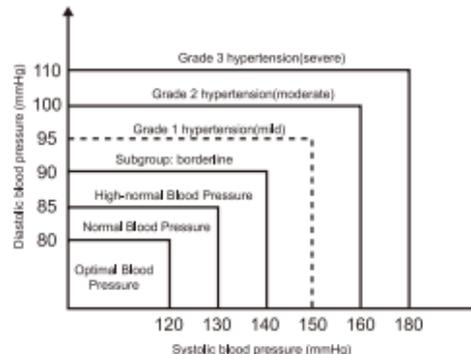

Blood Pressure (mm Hg) \ Level	Optimal	Normal	High-normal	Mild	Moderate	Severe
SYS	<120	120-129	130-139	140-159	160-179	≥180
DIA	<80	80-84	85-89	90-99	100-109	≥110

♥ **Irregular Heartbeat Detector**

An irregular heartbeat is detected when a heartbeat rhythm varies while the device is measuring systolic pressure and diastolic pressure. During each measurement, blood pressure monitor will keep a record of all the pulse intervals and calculate the average value of them. If there are two or more pulse intervals, the difference between each interval and the average is more than the average value of $\pm 25\%$, or there are four or more pulse intervals, the difference between each interval and the average is more than the average value of $\pm 15\%$, then the irregular heartbeat symbol will appear on the display with the measurement result.

CAUTION

The appearance of the IHB icon indicates that a pulse irregularity consistent with an irregular heartbeat was detected during measurement. Usually this is NOT a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that the device does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

♥ Why does my blood pressure fluctuate throughout the day?

1. Individual blood pressure varies multiple times everyday. It is also affected by the way you tie your cuff and your measurement position, so please take the measurement under the same conditions.
2. If the person takes medicine, the pressure will vary more.
3. Wait at least 3 minutes for another measurement.



♥ Why do I get a different blood pressure at home compared to the hospital?

The blood pressure is different even throughout the day due to weather, emotion, exercise etc. Also, there is the “white coat” effect, which means blood pressure usually increases in clinical settings.

What you need to pay attention to when you measure your blood pressure at home:

- If the cuff is tied properly.
- If the cuff is too tight or too loose.
- If the cuff is tied on the wrist.
- If you feel anxious.
- Taking 2-3 deep breaths before beginning will be better for measuring.
- Advice: Relax yourself for 4-5 minutes until you calm down.



♥ Is the result the same if measuring on the right wrist?

It is ok for both wrists, but there will be some different results for different people. We suggest you measure the same wrist every time.

This section includes a list of error messages and frequently asked questions for problems you may encounter with your blood pressure monitor. If the products not operating as you think it should, check here before arranging for servicing.

PROBLEM	SYMPTOM	CHECK THIS	REMEDY
No power	Display is dim or will not light up.	Batteries are exhausted.	Replace with new batteries
		Batteries are inserted incorrectly.	Insert the batteries correctly
Low batteries	+Lo Show on the display	Batteries are low.	Replace with new batteries
Error message	E 01 shows	The cuff is too tight or too loose.	Refasten the cuff and then measure again.
	E 02 shows	The monitor detected motion, while measuring.	Movement can affect the measurement. Relax for a moment and then measure again.
	E 03 shows	The measurement process does not detect the pulse signal.	Loosen the clothing on the wrist and then measure again.
	E 04 shows	Measure incorrectly.	Relax for a moment and then measure again.

TROUBLESHOOTING

PROBLEM	SYMPTOM	CHECK THIS	REMEDY
Error message	EExx, shows on the display.	A calibration error occurred. (XX can be some digital symbol, such as 01, 02, etc., if this similar situation appears, all belong to calibration error.)	Retake the measurement. If the problem persists, contact the retailer or our customer service department for further assistance. Refer to the warranty for contact information and return instructions.
Warning message	"out" shows	Out of measurement range	Relax for a moment. Refasten the cuff and then measure again. If the problem persists, contact your physician.

SPECIFICATIONS

Power supply	Battery powered mode: 2*AAA batteries (3V DC)
Display mode	Digital LCD V.A. 45mm x 33mm
Measurement mode	Oscillographic testing mode
Measurement range	Rated cuff pressure: 0mmHg ~ 299mmHg (0kPa ~ 39.9kPa) Measurement pressure: SYS: 60mmHg ~ 230mmHg (8.0kPa ~ 30.7kPa) DIA: 40mmHg ~ 130mmHg (5.3kPa ~ 17.3kPa) Pulse value: (40-199) beat/minute
Accuracy	Pressure: 5°C-40°C within ±3mmHg (0.4kPa) Pulse value: ±5%
Normal Working condition	A temperature range of: +5°C to +40°C A relative humidity range of 15% to 90%, non-condensing, but not requiring a water vapour partial pressure greater than 50 hPa An atmospheric pressure range of: 700 hPa to 1060 hPa
Storage & transportation condition	Temperature: -20°C to +60°C A relative humidity range of ≤ 93%, non-condensing, at a water vapour pressure up to 50hPa
Measurement perimeter of the wrist	About 13.5cm-21.5cm
Weight	Approx. 100g (Excluding the batteries)
External dimensions	Approx. 86mm × 66mm × 22mm (Excluding the cuff)
Attachment	2×AAA batteries, User manual, Quick start guide

SPECIFICATIONS

Mode of operation	Continuous operation
Degree of protection	Type BF applied part
Device Classification	Internally Powered ME Equipment
Protection against ingress of water	IP22: The first number 2: Protected against solid foreign objects of 12,5mm Ø and greater. The second number: Protected against vertically falling water drops when enclosure tilted up to 15°. Vertically falling drops shall have no harmful effects when the enclosure is tilted at any angle up to 15° on either side of the vertical.
Software Version	A03

WARNING: No modification of this equipment is allowed.

MANUFACTURER INFORMATION

♥ Manufacturer Information

Manufactured by: Harvard Medical Devices Ltd. HK

Company: Harvard Medical Devices Ltd. HK

Address: 1002, Railway Plaza, TST, HK

Authorized European Representative:

Company: Share Info Consultant Service LLC Repräsentanzbüro

Address: Heerdter Lohweg 83, 40549 Düsseldorf

♥ EMC Guidance

- 1) This product needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided, and this unit can be affected by portable and mobile RF communications equipment.
- 2) * Do not use a mobile phone or other devices that emit electromagnetic fields, near the unit. This may result in incorrect operation of the unit.
- 3) Caution: This unit has been thoroughly tested and inspected to assure proper performance and operation!
- 4) * Caution: This machine should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, this machine should be observed to verify normal operation in the configuration in which it will be used.

Guidance and manufacturer's declaration – electromagnetic emissions and Immunity

Table I

Guidance and manufacturer's declaration – electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Comply

Table 2

Guidance and manufacturer's declaration – electromagnetic immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for signal input/output 100 kHz repetition frequency	±2 kV for power supply lines Not Applicable 100 kHz repetition frequency
Surge IEC 61000-4-5	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV, ±2 kV common mode	±0.5 kV, ±1 kV differential mode Not Applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U_T ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T ; 1 cycle and 70% U_T ; 25/30 cycles; Single phase: at 0°. 0% U_T ; 250/300 cycle	0% U_T ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% U_T ; 1 cycle and 70% U_T ; 25/30 cycles; Single phase: at 0°. 0% U_T ; 250/300 cycle

Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC 61000-4-6	3 V 0, 15 MHz - 80 MHz 6 V in ISM and amateur radio bands between 0, 15 MHz and 80 MHz 80% AM at 1 kHz	3 V 0, 15 MHz - 80 MHz 6 V in ISM and amateur radio bands between 0, 15 MHz and 80 MHz 80% AM at 1 kHz
Radiated RF IEC 61000-4-3	10 V/m, 80 MHz - 2, 7 GHz 80 % AM at 1 kHz	10 V/m, 80 MHz - 2, 7 GHz 80 % AM at 1 kHz
NOTE U_T is the a.c. mains voltage prior to application of the test level.		

Table 3

Guidance and manufacturer's declaration - electromagnetic immunity								
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Maximum Power (W)	Distance (m)	IEC 60601-1-2 Test Level (V/m)	Compliance level (V/m)
	385	380-390	TETRA 400	Pulse modulation 18Hz	1.8	0.3	27	27
	450	430-470	GMRS 460, FRS 460	FM \pm 5kHz deviation 1kHz sine	2	0.3	28	28
	710	704-787	LTE Band 13, 17	Pulse modulation 217Hz	0.2	0.3	9	9
	745							
	780							
	810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18Hz	2	0.3	28	28
	870							
	930							

	1720	1700-1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217Hz	2	0.3	28	28
	1845							
	1970							
	2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28	28
	5240	5100-5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9	9
	5500							
	5785							